

Vet Temp 2

Veterinary Fluid Warmer

The easy way to prevent hypothermia!

- Use with standard fluids, blood and TPN fluids
- Light weight, 200 gr only!
- Portable
- Accepts standard tubing sets for small animals
- Fast and easy to set-up
- Easy to clean
- Over-temperature protection, five alarms including high and low temperature alarm
- Two attaching methods on I.V. stand lanyard and bracket



Vet Temp 2 supplies external heat to I.V. set fluids, near to the patient's body temperature.

The front and back of the device are joined by 2 hinges.

Simple construction and shock resistant casing for a high durability.

The operator can adjust the setting temperature in the range of 35°C - 42°C according to actual needs.



Vet Temp 2 Fluid Warmer



Why investing in a veterinary fluid warmer?

Maintaining the body temperature is very important during and after an anaesthesia for avoiding serious consequences such as:

- Risk of shock
- Risk of heart problems, renal insufficiency in a long term ...
- Risk of hypothermia

Vet Temp 2 helps to decrease the temperature difference between the animal and the I.V. liquid. The animal needs less energy for maintaining his body temperature.


Similar to the **Animat and Carbon Tech** heating Mats, the **Vet Temp 2** Fluid Warmer avoids the patient's temperature cools down too much.

Fields of applications and benefits:


- Safety infusion of a warm fluid
- Suitable for Infusions with a pump or for traditional infusions by gravity
- Easy to use
- Cost effective and durable device
- No need of disposables




Three lengths of heating groove design to make infusion heating more flexible and accurate:

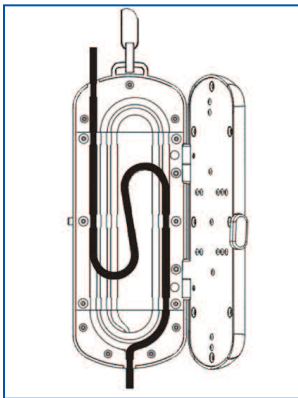

25 cm groove

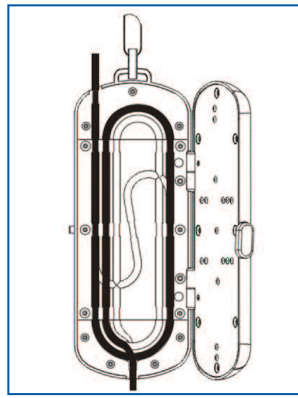

31 cm groove

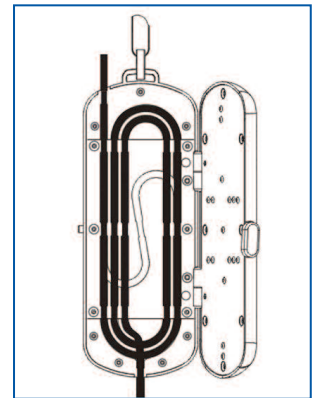

52 cm groove

Technical data:

Heating groove lengths:	25, 31 and 52 cm
Applicable I.V. set:	Exclusive IV set or standard PVC IV set (20 and 60 drops / ml)
I.V. set outer diameter:	Between 3.4 mm and 4.5 mm
Temperature:	Between 35°C and 42°C
Warming up time:	≤ 2 min.
Alarms:	Overheat & low temperatures, system error, door open, heating reminder
Power:	100-240 V / 50-60 Hz
Dimensions and weight:	18 x 7 x 3.5 cm - 200 gr



Medical equipment for veterinarians

15, rue du Noroit
22100 Taden - France
Tél. +33 (0)2 96 85 86 76
Fax +33 (0)2 96 85 86 77

welcome@manomedical.com

www.manomedical.com