



# Coag&Cut M120

**120 W Monopolar & Bipolar  
Electrosurgery unit**

**Perfect for daily small to medium surgeries**

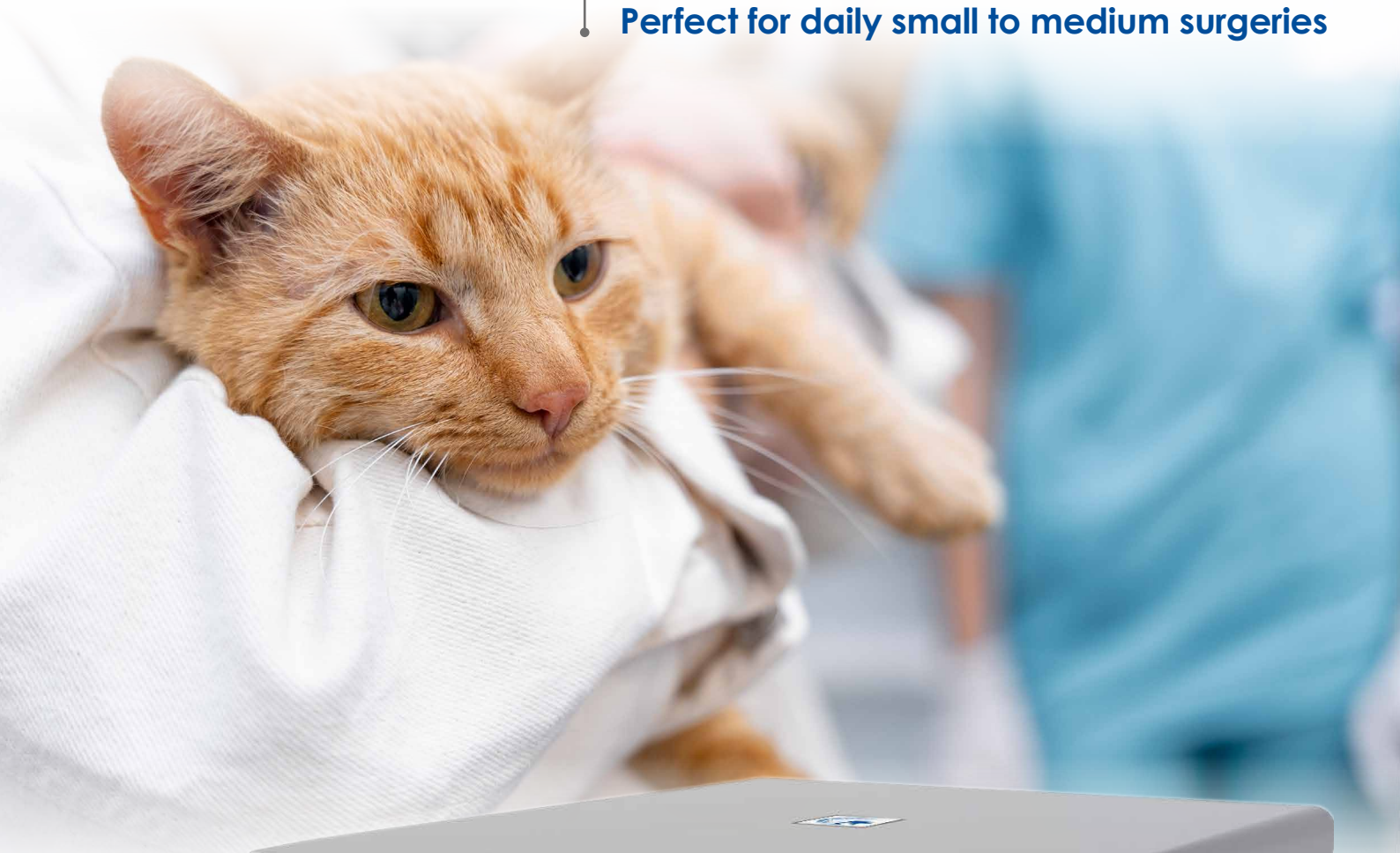

# Cut, coagulation and micro-coagulation

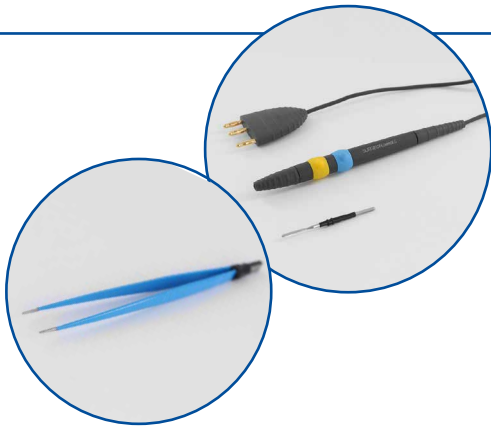

## • Perfect for daily small and medium surgeries

The high-performance **Coag&Cut M120** electrosurgery unit has a number of operating modes:

Cutting:

- ▶ CUT = pure Cut
- ▶ BLEND = Cut with Coagulation

Coagulation:

- ▶ FORCED = Superficial Coagulation
- ▶ SOFT = Deep Coagulation
- ▶ BIPOLEAR = Standard Bipolar Coagulation

## • Reliability

The digital power setting display and the continuous monitoring of the operational function by micro processor guarantee an absolute reliability of working conditions.

*Certified compliant MDD/93/42 EC Class: IIb*



## • Safety

The connection to the neutral plate is continuously monitored.

The monopolar handpiece, when used without a footswitch, allows to choose and activate either the Cut Mode (yellow button) or the Coagulation Mode (blue button). Coupled with the footswitch, it allows to choose the function which is activated only by the footswitch. Surgical procedures can be done safely.

*Compliant with safety standards EN60601-1-1 / EN60601-1-2 / EN60601-2-2*



## • Efficient

Maximum Output Power:

- ▶ CUT: 120W-250Ω
- ▶ BLEND: 90W-200Ω
- ▶ FORCED: 80W-150Ω
- ▶ SOFT: 60W-100Ω
- ▶ BIPOLEAR: 40W-100Ω

Working frequency: 600 kHz  
Main frequency: 50-60 Hz  
Electrical input power: 300 VA  
Size & Weight: 25 x 10 x 29 cm - 5 kg



## The Coag&Cut M120 package includes:

- ▶ A monopolar handpiece with a straight blade electrode
- ▶ A set of 10 additional electrodes (needles, loops, balls) + support
- ▶ A patient return electrode with cable
- ▶ A bipolar handpiece with cable
- ▶ A single footswitch

All the accessories are autoclavable to ensure durability and good hygiene



15, rue du Noroît  
22100 Taden - France  
Tel. +33 (0)2 96 85 86 76

[welcome@manomedical.com](mailto:welcome@manomedical.com)

[www.manomedical.com](http://www.manomedical.com)